EN: BREASTIUP – USER MANUAL

PRODUCT STANDARDS

Product Name: BREASTIUP

Usage Mode: A-F (6 modes)

Frequency: 1-100Hz

Operation Time: 12 minutes (per session)

Output: Up to 9.8mA maximum


INSTALLATION STEPS

1. Open the battery compartment.
2. Insert the batteries correctly (2 AAA batteries), then close the compartment.
3. Attach the abdominal muscle pad (or arm pad) to the main unit.
4. Remove the adhesive stickers from the pad and affix them separately onto the abdomen (or arms).

PRODUCT FEATURES

1. Low-frequency treatment involves rhythmic and gentle impulses. The specialized low-frequency electricity emits a regular rhythm or pauses. This rhythmical stimulation provides a soft, massage-like effect, akin to human hand pressure, kneading, and tapping.
2. **Six modes are available for selection:**
 - Easily choose from 6 modes.
 - The device automatically powers off after 12 minutes of inactivity to prevent energy wastage.

OPERATING INSTRUCTIONS

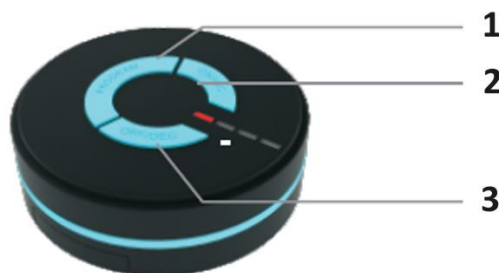1. ON: Power on / INC: Increase
2. OFF: Power off / DEC: Decrease
3. PROGRAM: Mode selection

When pressing the "power on" and "increase intensity" buttons, the indicator light gradually brightens. Conversely, pressing the "power off" and "decrease intensity" buttons dims the indicator light gradually.


- Intensity adjustment buttons, for both "power on" and "increase intensity," and "power off" and "decrease intensity," offer 10 levels of adjustment.

- Mode selection button: Choose from 6 modes (A-F).
- Battery compartment: Located at the rear for easy access.
- Dry cell batteries.

Mode Indicator Light

PAD USAGE INSTRUCTIONS

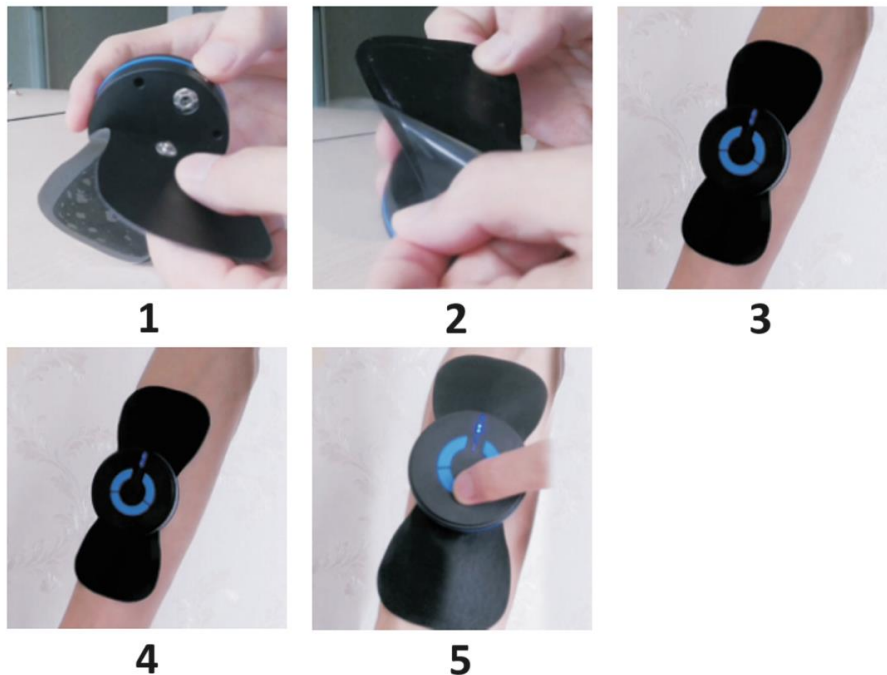1. Connect the main unit to the gel patch.
2. Remove the protective film from the gel patch.
3. Place the device on the desired massage area.
4. Press the "ON / INC" button to start the device; the light remains steady. Pressing the "ON / INC" button again makes the light flash.
5. Use the PROGRAM key to select the desired mode.
6. Adjust the massage intensity to your preference. Continue pressing the "ON / INC" button to increase intensity or press it to decrease intensity.

*Avoid direct contact with the adhesive patch immediately after activation.

MAINTENANCE

Main Unit:

Clean the unit's surface using a damp cloth and a neutral cleaning agent. Wipe dry with a clean cloth. Avoid storing the unit in high or low-temperature environments or areas with excessive humidity.

Adhesive Gel Pads:

If the adhesive side of the pad becomes dirty or loses its stickiness, it may cause discomfort. Use a damp towel to clean the dirty side for improved adhesion. Excessive moisture may affect adhesion. Replace adhesive gel pads (electrodes) with identical types and models.

BATTERY REPLACEMENT

1. Open the battery compartment.
2. Ensure correct insertion of batteries with regard to polarity (+ and -).
3. Close the battery compartment.
4. Turn on the switch and check the indicator light.

INSTRUCTIONS FOR RECYCLING AND DISPOSAL:

This label means that the product cannot be disposed of as other household waste throughout the EU. To prevent potential damage to the environment or human health from uncontrolled waste disposal. Recycle responsibly to promote the sustainable use of material resources. If you want to return a used device, use the drop-off and collection system, or contact the retailer from whom you bought the product. The retailer can accept the product for environmentally safe recycling.

A declaration by the manufacturer that the product complies with the requirements of the applicable EU Directives.

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