## EN: CUPPY – USER MANUAL

### BASIC PARAMETERS

Charging voltage DC 5V

Working voltage: 3.7V

Power: 8W

### PRODUCT PARTS

A – Power button

B – Heating button

C – Negative pressure button

D – LCD screen

E – Time button

### FEATURES

#### Negative pressure gravitation

- A. Aerobic exercise for sucking up the oil tissue of the skin.
- B. Increase lymphatic function and blood circulation.
- C. Excretion of toxins through the lymphatic circulation.
- D. Scraping effect.

#### Thermal energy

- A. The cortical collagen is grown by the thermal effect.
- B. Continuous heat, relax the spirit, release the pressure.

#### Magnetic wave therapy

- A. Help the swelling subsided, eliminate swelling and pain.
- B. Eliminate fatigue, help physical recovery.

### MASSAGE TECHNIQUES

1. Neck: Brush from top to bottom, often wipe the neck.
2. Back: Scrape from top to bottom.
3. Lower limbs: Dredge from top to bottom.
4. Foot sole: Massage from the heel to the toe.
5. Abdomen: The abdomen is dredged from top to bottom and can also be operated by rotating the circular motion.
6. Upper limbs: Dredge from top to bottom.

7. Legs: Scrape from the lower leg to clear the lymphatic toxins in the legs and shape the legs. Eliminates edema, eliminates dry legs, and eliminates the choice of leg fat.

## PRECAUTIONS

- If your skin has various veins, it should not be used.
- Do not use if pregnant or during period.
- Do not use directly after a meal.
- Do not expose the treated skin part to any exposures such as water, lotion, soaps, sun,..
- When using, follow the same direction.
- Use only with daily intervals, depending on skin condition.

Note: The device has a hot surface. Care must be taken when using.


## HOW TO USE

1. After cleansing the skin, apply a small amount of essential oil.
2. Fully charge the instrument first.
3. Press the power button and the digital light will be on. At this time, you can adjust the "heating" and "negative pressure" according to the customer's needs.
4. Press the "heat" button, the temperature is divided into 1-9 settings, each up, the higher the temperature.
5. Press the "negative pressure" button, "negative pressure" is divided into 1-9 settings. For each upward adjustment, the negative pressure is higher.
6. After using the instrument, press the middle power button to turn off the instrument, clean the instrument, and place it in a dry place.

Note: The default working time of the instrument is 20 minutes. After the time passes, the instrument will stop working automatically for 20 minutes. This will extend the lifespan of the product.

## IMPORTANT NOTE

Before using the instrument, you must put the filter cotton into the filter plug (Figure 1), and then put the filter plug into the suction hole of the instrument (Figure 2). It can effectively block the essential oil or massage products directly into the instrument. Effectively extend the life of the instrument. (Replace the filter cotton after each use, the filter cotton is air-dried and can be reused).
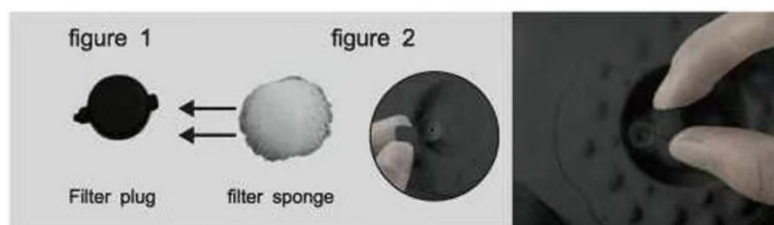
## INSTRUCTIONS FOR RECYCLING AND DISPOSAL:

This label means that the product cannot be disposed of as other household waste throughout the EU. To prevent potential damage to the environment or human health from uncontrolled waste disposal. Recycle responsibly to promote the sustainable use of material resources. If you want to return a used device, use the drop-off and collection system, or contact the retailer from whom you bought the product. The retailer can accept the product for environmentally safe recycling.

A declaration by the manufacturer that the product complies with the requirements of the applicable EU Directives.

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