


EN: GROWNETIC – USER MANUAL

PRODUCT OVERVIEW

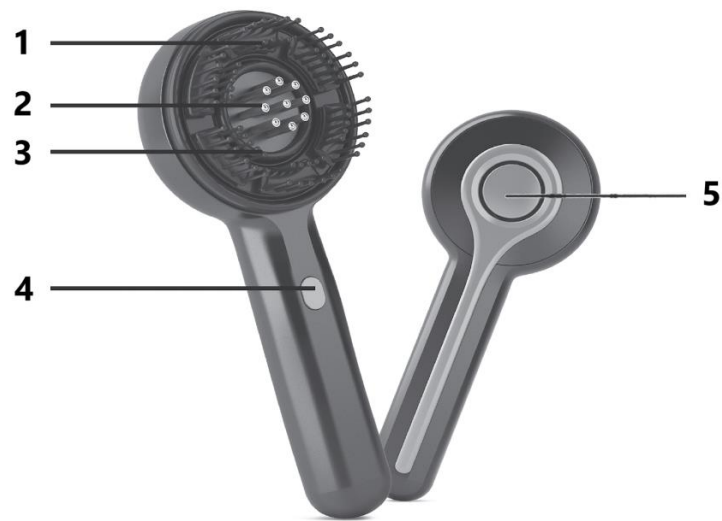


1. Massage Head
2. Fluid Guide Head
3. Infrared Light
4. On/Off Button
5. Liquid Tank Cover

Product Specifications

- **Product Name:** GROWNETIC
- **Rated Voltage:** 3.7V
- **Rated Power:** 5W
- **Battery Capacity:** 1400mAh
- **Charging Interface:** Type-C
- **Package Size:** 242 x 116 x 65 mm

Product Functions

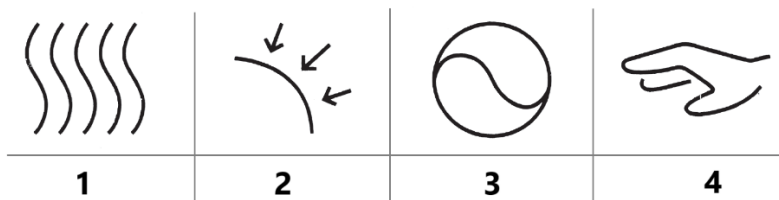


1. Scalp Massage
2. Infrared Light Therapy
3. Ball-Draining Liquid Dispersion
4. Ergonomic Grip

Packaging Contents




1

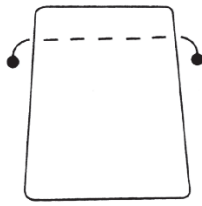

2



3

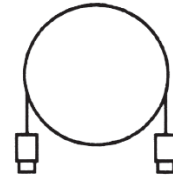

4

1. Main Device
2. Storage Bag
3. User Manual
4. Charging Cable

SAFETY PRECAUTIONS

1. This product is a household electronic healthcare device and cannot replace medical equipment. It is not intended for diagnosis, treatment, or other commercial purposes.
2. Remove accessories or objects that may interfere with the massage before use.
3. Discontinue use if you experience pain or discomfort in the treated area.
4. If any irregularities occur while using the device, turn it off, unplug it, and contact customer support. Do not attempt self-repair.
5. Follow all safety guidelines when using the product in specific environments.
6. Use only compatible charging equipment with the correct voltage to avoid damage. Do not mix with other charging devices.
7. Keep the product away from fire, moisture, and corrosive environments.

HOW TO USE

Note: This product cannot be used while charging. Always unplug the power adapter before use.

Powering On & Mode Selection

- Long press the On/Off button for about **2 seconds** until a "**beep**" sound is heard, indicating power on. The device starts in **Mode 1** by default.
- Press the On/Off button to switch between modes. A "**drip**" sound confirms a successful mode change.
 - **Mode 1:** Fast Massage + Infrared Light
 - **Mode 2:** Slow Massage + Infrared Light
 - **Mode 3:** Alternating Fast & Slow Massage + Infrared Light

Adding Liquid

1. Open the liquid tank cover and fill it with a suitable hair care product. Do not exceed **3ml** to avoid waste.
2. Clean the liquid storage bin after each use to maintain hygiene.

3. Different care solutions may affect the discharge rate. Lower concentrations may cause slight leakage—clean promptly if needed.
4. **Warning:** Do not add unknown or corrosive substances to prevent damage.

Smart Timing Function

The device has a **10-minute** automatic shutoff function.

Powering Off

Long press the On/Off button to turn off the device.

Cleaning Instructions

1. To clean the outer surface, rinse it with water. **Do not expose the charging port to water.**
2. Rinse the liquid storage bin and outlet hole with clean water after each use to prevent blockages.
3. Wipe the device dry and store it in a **ventilated and dry place.**

OTHER INFORMATION

1. Normal operation may produce sounds such as "hum" or "click," which are not malfunctions.
2. New products may contain residual liquid in the tank due to factory quality testing. This is normal and safe.

CHARGING INSTRUCTIONS

1. Connect the **Type-C** charging cable to the charging port.
2. Plug the other end into a USB charging port or a **5V adapter** and connect it to a power socket.
3. During charging, the indicator light flashes. When fully charged, the red light remains on.
4. Charging time is approximately **3 hours.**

Notes:

- Ensure the charging port is dry before charging.
- Slight heating of the device and adapter during charging is normal.

End of Charging

Unplug the adapter and remove the charging cable.

TROUBLESHOOTING

If you experience an issue, check the following before contacting support:

Issue	Possible Cause	Solution
Device won't turn on	Low battery	Charge the device
Weak massage strength	Low battery	Charge the device


Important Notices

- Before disposal, remove the battery following regulations.
- Always **turn off** the device before removing the battery.
- Dispose of the battery in compliance with local environmental laws.
- **Non-professionals must not attempt disassembly.**

BATTERY REMOVAL PROCEDURE

1. Disconnect the external power source and remove the device casing.
2. Unplug and carefully remove the battery.
3. Dispose of the battery at a designated recycling facility.
4. Follow local environmental regulations for proper disposal.

OPENING THE COVER

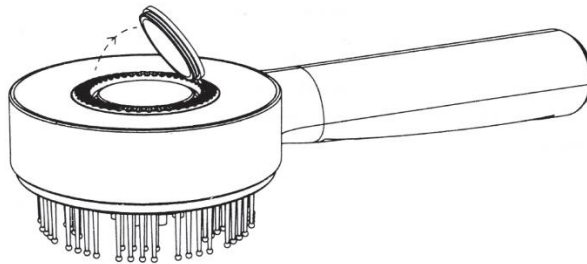

STORAGE & TRANSPORTATION

- If storing for **more than 3 months**, recharge the device before use.
- Do not place the device near walls, floors, or ceilings to prevent damage.
- Ensure the storage area is **ventilated, dry, and free from direct sunlight or corrosive gases.**
- The product is designed with shock protection suitable for **air, rail, road, and sea transport.**
- Avoid exposure to **rain, snow, inversion, or impact** during transportation.

POWER CHECK & MAINTENANCE

- Always verify that the power supply voltage is within the specified range before use.
- Protect the device from moisture, dust, and rodents.
- Clean the device **only when powered off.**
- Store the device in a **safe place** and avoid placing heavy objects on it.
- In case of malfunction, seek professional repair services. **Do not attempt self-repair.**

INSTRUCTIONS FOR RECYCLING AND DISPOSAL:



This label means that the product cannot be disposed of as other household waste throughout the EU. To prevent potential damage to the environment or human health from uncontrolled waste disposal. Recycle responsibly to promote the sustainable use of material resources. If you want to return a used device, use the drop-off and collection system, or contact the retailer from whom you bought the product. The retailer can accept the product for environmentally safe recycling.



A declaration by the manufacturer that the product complies with the requirements of the applicable EU Directives.