
PRODUCT INTRODUCTION

Parts overview

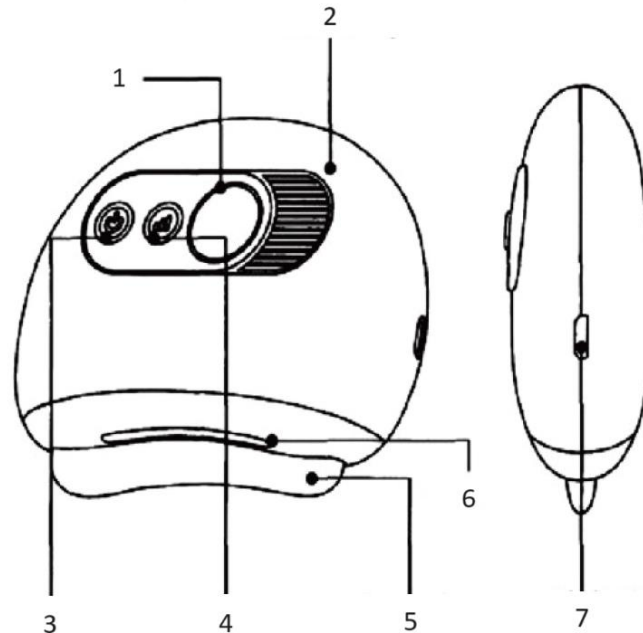1. Display screen
2. Machine casing
3. Switch / Mode button
4. Gear adjustment button
5. Scraping instrument board
6. Red light
7. USB charging port


OPERATION METHOD

Press and hold the switch/Mode button for 2 seconds to enter the power on state. The indicator board is on red Starting, temperature default 1, vibration default 0. Press the switch/mode key, select the function, temperature and vibration switch in turn, at the corresponding function position, the corresponding mode indicator will light up (temperature mode function indicator is red, vibration mode function indicator is blue), adjust the gear through the gear key, vibration 0-9, temperature gear 0-9. (1-3 is low grade, 4-6 is medium grade, 7-9 is high grade) Long press the gear key for 2 seconds to turn the red light off/on. Press and hold the switch/Mode button for 2 seconds to shut down. Low voltage warning: The display blinks for 5 times and shuts down automatically. Display 8 flickers when charging, while display 8 lights up when full. Automatic shutdown after 30 minutes of work.

Note: The product cannot be turned on when charging, and can only be turned on by unplugging the power supply. Warning: The appliance has a hot surface, and people who are not sensitive to heat must pay attention when using it,

PRODUCT USAGE DIAGRAM

In general, the order of use is head, neck, back (thoracic vertebrae, Lumbar vertebrae, sacra vertebrae), thoracic abdomen, upper limbs, lower limbs.

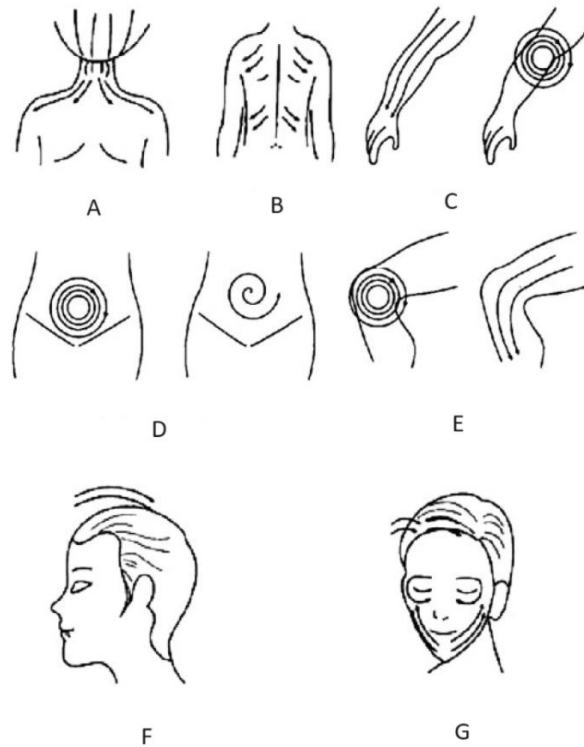- A. Shoulder massage
- B. Back massage
- C. Elbow massage
- D. Abdominal massage
- E. Leg massage
- F. Head massage
- G. Facial massage

CONTRAINDICATION

1. Critical conditions, such as severe heart disease, kidney failure, cirrhosis of the liver, ascites, complete do not use for severe edema.
2. People with highly sensitive or severe skin allergies are prohibited from contact skin diseases.
3. Do not use this product where there are lumps, ulcerations, macules and unexplained masses on the body surface
4. People with bleeding tendencies, such as late diabetes, severe anemia, leukemia, a plastic anemia and thrombolytic are strictly prohibited.

5. Low blood pressure, low blood sugar, excessive weakness and nervous tension, especially fear of pain Do not use.
6. People who have been ill for a long time, are extremely weak, thin, and the skin has lost its elasticity are strictly prohibited
7. It is strictly prohibited to use the painful parts of acute sprains, wounds or fractures.
8. Over-thirst, over-hunger, over-drinking, drunkenness, over-fatigue, and over-staying up late are severe Do not use.
9. Do not use the lower abdomen during menstruation; Prohibited for pregnant women.
10. It is prohibited for patients with skin perception disorders or those who are not sensitive to heat.
11. This product should not be used by postoperative patients and those who have been fitted with a pacemaker.

PRODUCT PRECAUTIONS AND SAFETY

1. Pay attention to the stone needle temperature when using, choose your own appropriate temperature gear,The temperature acceptable to the skin shall prevail. Do not overheat to prevent burns.
2. When the smell remains after use with essential oil, wipe the product with a clean wet cloth to keep the product clean. Before giving it to others, 75% alcohol should be used for disinfection to avoid infection.
3. Do not move the instrument by holding the power cord.
4. This product can not fall from a high altitude, impact, should be handled gently, careful Fragile.
5. This machine is not waterproof, do not soak in water or water.
6. It is strictly prohibited to use chemical solvents to clean the body, power supply and key cover of the machine Board.

Note: This product cannot be used as a substitute for drugs or medical devices.

PRODUCT PARAMETERS

Product name: GUASHAGLOW

Product model: LSF-888

Rated input voltage: 5V

Rated input power: 5W

RECHARGEABLE BATTERY REMOVAL METHOD

- The battery must be removed from the appliance before it is discarded;
- The appliance must be powered off when the battery is removed;
- The battery should be disposed of safely.

1. The built-in rechargeable lithium battery of this product may pollute the environment When discarding this product, please do not take it with general life The garbage should be put together and the battery should be taken out.

2. Before decomposing the product and removing the lithium battery, be sure to turn off the power or power off.

3. Take out the steps of the lithium battery, pry open the front and back covers of the main body, so that the front and back covers are separated, cut the battery lead with scissors, and remove the battery.

4. The removed lithium batteries should be disposed of safely and put into the battery recycling center designated by the government, which is conducive to environmental protection.

INSTRUCTIONS FOR RECYCLING AND DISPOSAL:


This label means that the product cannot be disposed of as other household waste throughout the EU. To prevent potential damage to the environment or human health from uncontrolled waste disposal. Recycle responsibly to promote the sustainable use of material resources. If you want to return a used device, use the drop-off and collection system, or contact the retailer from whom you bought the product. The retailer can accept the product for environmentally safe recycling.

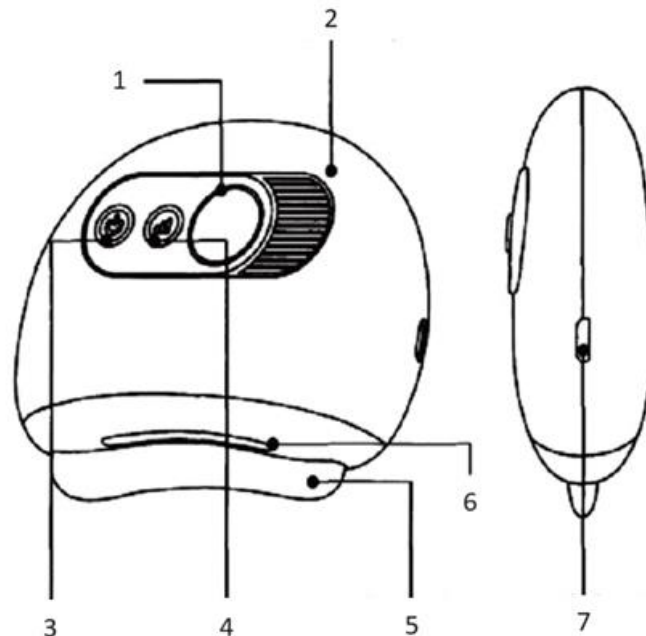A declaration by the manufacturer that the product complies with the requirements of the applicable EU Directives.

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