

EN: LUMEVAGLOW – USER MANUAL

PRODUCT PARAMETERS

Product name: LUMEVAGLOW

Product model No.: ES-1081

Charging voltage: 5V, 1A

Rated power: 2.5W

Battery model No.: 602535

Battery parameters: 3.7V, 500mAh

Charging time: 2 hours

Product size: 89x47x120mm

Executive standard: GB 4706.1-2005, GB 4706.15-2008

Working time: The device enters standby mode after 4 minutes and automatically shuts off after 1 minute of inactivity.

FUNCTION INTRODUCTION

This product is designed for skin cleansing, EMS care, and moisturizing at a warm temperature. It integrates positive and negative ions, high-frequency vibration, red light, blue light, and green light.

- **Positive and Negative Ions:** Using the principle of opposite charge attraction, bioelectric currents help extract deep-seated dirt from pores, while cosmetic ingredients penetrate deep into the skin.
- **Vibration:** Enhances skin tightness, elasticity, and metabolism through mechanical vibration and voltage stimulation.
- **Warm Massage:** Promotes absorption, improves skin texture, and gives a plump, youthful look.
- **Red Light:** Penetrates deeply to revitalize collagen, brighten the skin, and reduce fine lines and wrinkles.
- **Blue Light:** Stimulates collagen synthesis, rejuvenates the skin, tightens sagging areas, refines pores, and improves oily, sensitive, and acne-prone skin.
- **Green Light:** Enhances oxygen utilization, promotes skin microcirculation, accelerates lymphatic drainage, reduces puffiness, and soothes the skin.

PRODUCT CHARACTERISTICS

CLEAN Mode (Cleaning and Exporting): Removes deep-seated dirt, shrinks pores, fades discoloration, and helps treat acne.

EMS Mode (EMS Nursing): Enhances skin absorption of skincare products, boosts skin activity, and supports skin rejuvenation, whitening, and wrinkle reduction.

HOT Mode (Warm Import): Stimulates circulation, improves skin tightness and elasticity, and maintains skin hydration for a youthful glow.

PRODUCT COMPONENTS



1. Rejuvenation device head
2. Red, blue, and green light
3. Gear indicator light
4. Function indicator light
5. Function regulation (On/Off)
6. USB charging interface

NOTES FOR USE

1. Long press the power button for 1.5 seconds to turn on the device. It cannot be used while charging.
2. Hold the induction bars on both sides during use to complete the circuit with the body.
3. The device uses microcurrent ion technology, so a light tingling sensation is normal.
4. When extracting deep-seated dirt, use makeup remover or cleansing water—do not use facial cleansers or cleansing oils.
 - First, use the device to draw out deep-seated impurities, then wash off surface dirt with a facial cleanser.

SAFETY PRECAUTIONS

1. Not recommended for individuals who have undergone facial plastic surgery, have metal implants, or are pregnant.
2. Wipe the introduction head with a tissue. Do not wash with water or touch metal surfaces.
3. Keep the device away from bathrooms and humid environments to avoid short circuits.

4. Use each mode only once per day, with a recommended three-mode sequence per use.
5. The device is designed for facial and neck use only—avoid areas with acne, blemishes, or redness.
 - Those with sensitive skin should use it only once a week.
6. Ensure the power is off before performing maintenance. Cleaning with a wet towel while powered on may cause mild electric shocks.
7. Do not use corrosive solutions (e.g., alcohol, acetone, gasoline) to clean the device, as they may cause malfunctions, cracks, or discoloration.

USAGE INSTRUCTIONS

CLEAN Mode (Cleaning and Exporting)

- Effectively removes pore dirt and enhances skin translucency. Avoid thick cotton pads.
- Combined with green light to promote oxygen utilization, improve microcirculation, and soothe the skin.

Steps:

1. Hold the power button for 1 second to enter "CLEAN" mode (green light indicator).
2. Apply makeup remover or cleansing lotion evenly onto the skin. (Device is not waterproof—do not spray water directly on it.)
3. Hold the metal induction bars on both sides to complete the ion circuit.
4. Glide the device upward along the face and neck.

EMS Mode (Electromuscular Stimulation)

- Stimulates and relaxes the skin, helping skincare products penetrate deeply.
- Combined with red light for collagen revitalization, skin brightening, and wrinkle reduction.

Steps:

1. Hold the power button for 1 second to start, then press again to enter "EMS" mode (red light indicator).
2. Apply toner or lotion evenly onto the skin.
3. Hold the metal induction bars to complete the ion circuit.
4. Glide the device upward along the face and neck.


HOT Mode (Warm Import)

- Gently wakes and relaxes the skin, restoring radiance and elasticity.
- Combined with blue light to stimulate collagen synthesis, refine pores, and treat oily, sensitive, or acne-prone skin.

Steps:

1. Hold the power button for 1 second to start, then press twice to enter "HOT" mode (blue light indicator).
2. Apply a facial mask or serum evenly onto the skin.
3. Hold the metal induction bars to complete the ion circuit.
4. Glide the device upward along the face and neck.

RECOMMENDED MOVEMENTS

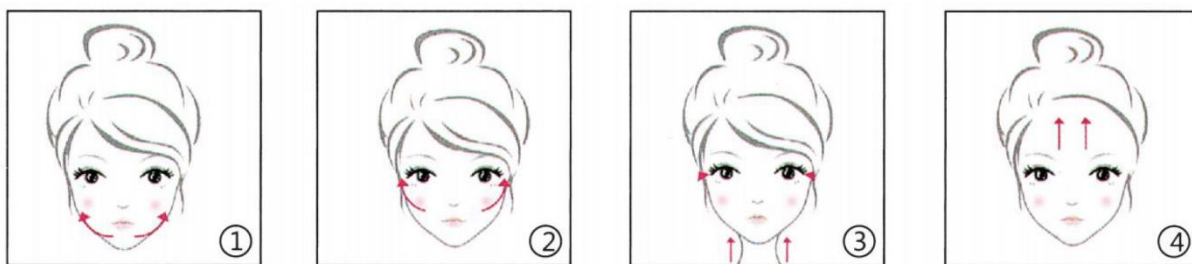

1. Lift from chin to cheek.
2. Lift from cheek to outer eye area.
3. Lift from neck to chin.
4. Lift from chin to jawline.

INSTRUCTIONS FOR RECYCLING AND DISPOSAL:



This label means that the product cannot be disposed of as other household waste throughout the EU. To prevent potential damage to the environment or human health from uncontrolled waste disposal. Recycle responsibly to promote the sustainable use of material resources. If you want to return a used device, use the drop-off and collection system, or contact the retailer from whom you bought the product. The retailer can accept the product for environmentally safe recycling.



A declaration by the manufacturer that the product complies with the requirements of the applicable EU Directives.