EN: MASERKO – USER MANUAL

Dear Customer,


Thank you for purchasing the MASERKO massage mat. This product has been rigorously tested to ensure safety and reliability. Please read this manual carefully before using the product to understand its proper use. Pay special attention to the safety warnings, and keep this manual for future reference.

Note: Our company reserves the right to change the product design and details without prior notice. The actual product may vary from the illustrations provided in this manual.

PRODUCT PARAMETERS


- **Name:** MASERKO
- **Voltage:** 12V DC
- **Model:** HS-D5
- **Power:** 24W
- **Standards:** GB4706.1-2005, GB4706.10-2008

SCHEMATIC DIAGRAM OF PRODUCT

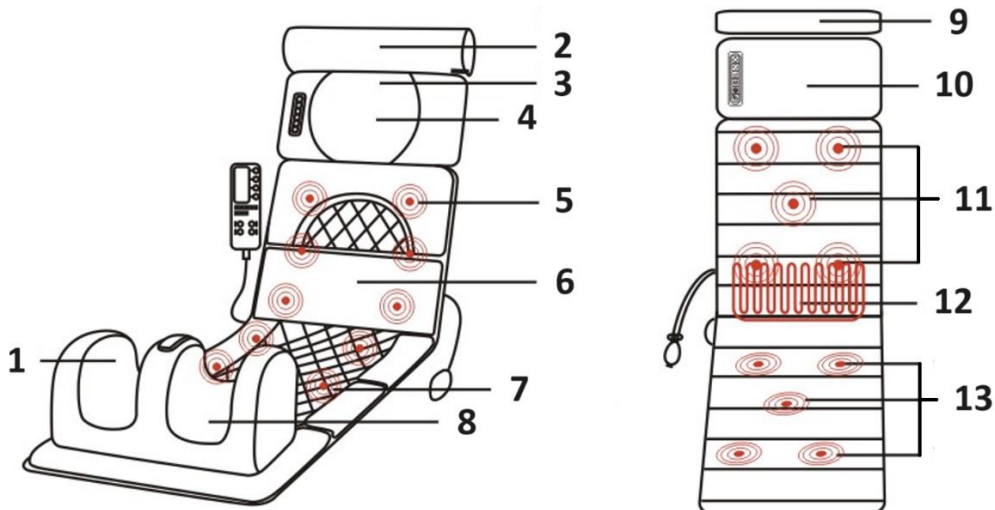1. Leg Massager
2. Multiple Herb Pillow
3. Massage Pillow
4. Head Massager
5. Vibration Point
6. Airbag
7. Vibration Point
8. Leg Massager
9. Tsao Stick (Moxibustion)
10. Head Massage
11. Vibration Points
12. Heater
13. Vibration Points

MANUAL CONTROLS OVERVIEW

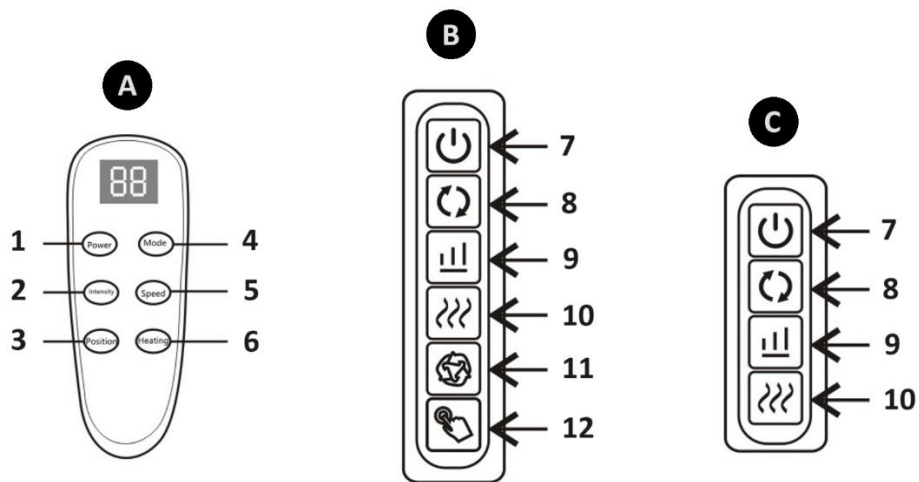A – Corded Controller

B – Head Massager Control Buttons (on pillow)

C – Leg Massager Control Buttons (on leg rest)

1. Power Button
2. Intensity Switch
3. Position Control
4. Mode Switch
5. Speed Adjustment Switch
6. Heating Adjustment Switch
7. On/Off Button
8. Turn/Inversion
9. Intensity Switch
10. Hot Compress
11. Energy-saving Mode
12. Smart Mode

USER INSTRUCTIONS

- **Switch/Timing:** Press the Power button on the controller to switch on. Press again to set the timer to 10, 20, or 30 minutes.
- **Mode Selection:** Press the Mode button on the controller to select from 5 different massage modes.
- **Intensity Adjustment:** Press the Intensity button on the controller to choose from 3 levels of massage intensity.
- **Position Adjustment:** Press the Position button on the controller to target the shoulders, back, arms, or legs for a customized massage.
- **Waist Airbag:** For models with a waist airbag, press the balloon to inflate it for additional comfort. Deflate when not in use.

BASIC SAFETY INSTRUCTIONS

- Disconnect the plug after use or during cleaning.
- Do not use in bathrooms or wet areas. Store in a dry place away from direct sunlight.
- Do not use while sleeping to avoid prolonged stimulation of one area.
- Do not use the massage mat within half an hour after a meal to prevent adverse effects.
- Stop use immediately if you experience discomfort due to excessive stimulation.
- Children or individuals with physical limitations should use the product under supervision.
- Consult a doctor before use if you are pregnant, have spinal pain, acute diseases, heart conditions, or other medical issues.
- Drink 500cc of boiled water or hot tea after massage to help expel metabolites from the body.
- If the power cord is damaged, contact the manufacturer for a replacement to avoid safety hazards.
- If the surface becomes hot during use, those with reduced heat sensitivity should be cautious.
- This product is not intended for use by persons with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless supervised.
- Supervise children to ensure they do not play with the product.

BASIC FAULT ANALYSIS / TROUBLESHOOTING

Common Fault Symptoms and Possible Causes

- **The Product Stopped Working**
 1. Is the adapter plugged into the power outlet?
 2. Is the hand controller turned on?
 3. Is the DC plug of the adapter connected to the DC plug of the massage mat?

Handling Method:

1. Ensure the adapter is firmly plugged into the power outlet.
2. Turn on the hand controller.
3. Connect the DC plug of the adapter to the DC plug of the massage mat.

- **The Product is Noisy**

Is it being used on a hard surface?

Handling Method:

Place the product on a soft bed, sofa, or use a cushion underneath.

- **Low Intensity**
 1. The grid voltage is too low.
 2. The massage area is not targeted properly.

Handling Method:

1. Use a voltage stabilizer or wait until the voltage stabilizes.
2. Ensure the massage part is accurately targeted.

If these methods do not resolve the issue, please contact your local dealer or our company. If the massage mat malfunctions or cannot be used, unplug the adapter. Users should not disassemble or repair the massage mat themselves, as this will void the warranty and our company will not be responsible for any resulting damage.

INSTRUCTIONS FOR RECYCLING AND DISPOSAL:

This label means that the product cannot be disposed of as other household waste throughout the EU. To prevent potential damage to the environment or human health from uncontrolled waste disposal. Recycle responsibly to promote the sustainable use of material resources. If you want to return a used device, use the drop-off and collection system, or contact the retailer from whom you bought the product. The retailer can accept the product for environmentally safe recycling.

A declaration by the manufacturer that the product complies with the requirements of the applicable EU Directives.

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