
EN: PEDIPRO – USER MANUAL

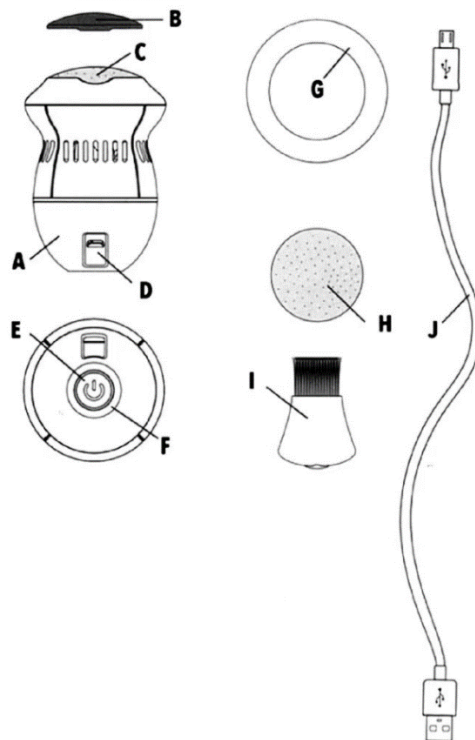- A- PEDIPRO (main unit)
- B- Emery Pad Cover
- C- Emery Pad
- D- Charging Port
- E- Power Button & Low/High Selection
- F- Battery Indicator & Speed Light
- G- Circular Sponge Filter
- H- Buffing Pad
- I- Emery Pad Cleaning Brush
- J- USB Charging Cable

HOW TO USE

1. Start by removing the emery pad cover by lifting up on one side.
2. Plug the micro-USB connector end of the supplied USB cable into the charging port of the PEDIPRO. Plug the other end of the USB cable into a USB wall adapter, then plug the adapter into the wall outlet. Turn on PEDIPRO by holding down the power button for 2 seconds. To turn off PEDIPRO hold down the power button until light turns off.
3. When unit is in the on-position, speed is low. To achieve high speed, click the power button once.
4. Holding the PEDIPRO securely, gently press the emery pad to the thickest area of dead skin you want to remove and slowly move outwards towards healthy skin.
5. Check intermittently for desired results. Once new, healthy skin is revealed or callus is removed, stop use and switch the PEDIPRO off.

HOW TO CLEAN

1. Start by removing the emery pad. Emery pad can be cleaned by running under warm water. Allow the emery pad to air dry. Do not attempt to clean the main unit by running under water. To clean the main unit wipe down with a soft damp cloth.
2. Once emery pad is removed, turn the PEDIPRO cover counterclockwise to open shavings chamber.
3. Slowly remove circular sponge filter from chamber. Rinse filter with warm water and allow to air dry. NOTE: Circular Sponge Filter and Emery Pad can be used up to 100x until replacement is needed.
4. Clean the remaining skin shavings out of the chamber with the provided shavings remover brush.
5. Once chamber is fully cleaned replace the circular filter, shavings chamber cover & emery pad securely. NOTE For best performance, clean PEDIPRO after each use.

WARNING: Do not immerse the unit in water or allow water anywhere near the unit. PEDIPRO is NOT waterproof. Store in a dry place. Such measures reduce the risk of electric shock and damage.

INSTRUCTIONS FOR RECYCLING AND DISPOSAL:


This label means that the product cannot be disposed of as other household waste throughout the EU. To prevent potential damage to the environment or human health from uncontrolled waste disposal. Recycle responsibly to promote the sustainable use of material resources. If you want to return a used device, use the drop-off and collection system, or contact the retailer from whom you bought the product. The retailer can accept the product for environmentally safe recycling.

A declaration by the manufacturer that the product complies with the requirements of the applicable EU Directives.

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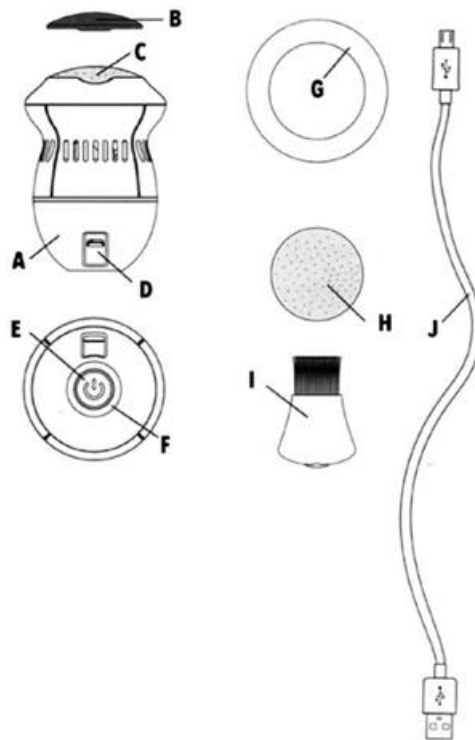- A- PEDIPRO(glavna enota)
- B-Pokrovček za šmirgljevo blazinico
- C-Smirkova blazinica
- D-Polnilna vrata
- E-Gumb za vklop in izbira nizke/visoke vrednosti
- F-Indikator baterije in lučka hitrosti
- G-Krožni gobast filter
- H-Blazinica za brušenje
- I-Čopič za čiščenje s smirkovo blazinico
- J-Polnilni kabel USB

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