EN: SIX PACK – USER MANUAL

Thank you for choosing the SIX PACK muscle stimulator. This manual has been sent to your email along with the invoice, and it's also available on our website under the 'User Manual' tab.

BASIC SAFETY INSTRUCTIONS

Before use, please carefully read and keep this manual. It contains important safety information. The instructions are divided into three categories to guide you in handling and using the device safely.

Individuals with the following conditions should not use this product:

- Those using electric medical equipment
- Individuals who have undergone plastic surgery on the targeted body parts
- Pregnant women
- Individuals who have recently undergone surgery

Incorrect operation of the device can pose risks.

WHAT IS LOW-FREQUENCY TREATMENT?

The body naturally generates a low level of electricity, which is imperceptible to us. However, external electronic impulses can stimulate sensitivity, leading to various bodily responses. This stimulation can help alleviate fatigue, reduce stress, and provide other benefits, such as:

- Alleviating tightness in the shoulders, waist, and back
- Improving muscle paralysis
- Providing a massage-like effect
- Aiding in weight management


PRODUCT SPECIFICATIONS

- Product name: Six Pack
- Battery: 1.5V AAA x 2 (not included)
- Modes: 6 (A-F)
- Frequency: 1 to 100Hz
- Operating time: 12 minutes per session
- Output: up to 9.8mA
- Dimensions: 21 x 10 x 2.3cm
- Weight: 57.4g

PRODUCT FEATURES

1. The low-frequency treatment delivers gentle and rhythmic impulses. These impulses mimic the effects of a massage, offering comfort and relaxation.
2. Choose from 6 modes for personalized treatment (A-F).


MAIN UNIT OVERVIEW

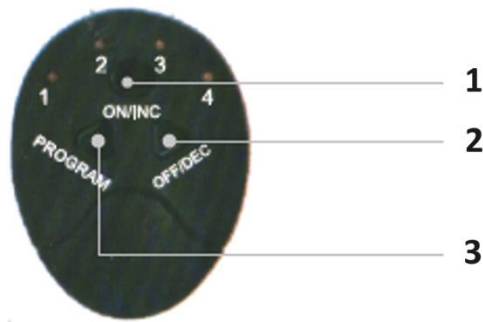1. ON: Power on / INC: Increase intensity
2. OFF: Power Off / DEC: Decrease intensity
3. PROGRAM: Mode switch

- Press ON to power on/increase intensity, and OFF to power off/decrease intensity.
- Use PROGRAM to switch modes.
- The device automatically powers off after 12 minutes to prevent overuse.

Mode indicator

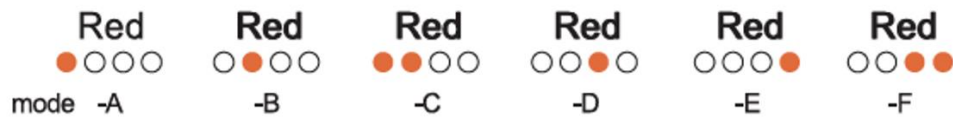MAIN UNIT MAINTENANCE

1. Clean the unit's surface with a damp cloth using a mild detergent, then wipe dry.
2. Store the device in a moderate temperature and humidity environment.

HOW TO USE THE PADS

1. Clean the intended application area before use.
2. Remove the adhesive gel pads from the vinyl cover and apply them to desired body parts.
3. Ensure the device is turned off before attaching the pads.
4. Avoid touching the adhesive side to prolong usage.
5. Replace pads as needed with the same type and model.

PRODUCT OPERATION

1. Begin with the lowest intensity setting.
2. Apply the gel pads to desired areas.
3. Power on the device.
4. Select a mode.
5. Adjust intensity as needed.
6. Each session lasts 12 minutes.

7. Power off using the OFF/DEC button when finished.
- Limit sessions to 12 minutes, 2-3 times per day.
 - Avoid use immediately after meals; preferable after a shower or warm wash.

POST-USE CARE

1. Always power off after use.
2. Remove gel pads from the machine and attach them to the pad holder.
3. Clean the device with a damp cloth and mild detergent, then wipe dry.
4. Store in a cool, dry place away from direct sunlight, heat, humidity, and dust.
5. Keep out of reach of children.
6. Avoid touching battery terminals with wet hands.
7. If not in use for more than a month, remove batteries and store separately.

CHANGING THE BATTERY

1. Open the battery cover.
2. Insert two AAA batteries correctly.
3. Close the cover securely.
4. Turn on the device and ensure the indicator light is functioning.

INSTALLATION INSTRUCTIONS

1. Open the battery compartment.
2. Insert two AAA batteries correctly.
3. Attach the device to the abdominal pad.
4. Peel off the adhesive pads.
5. Attach the device to the arm pads.
6. Apply the pads to your arms.

INSTRUCTIONS FOR RECYCLING AND DISPOSAL:

This label means that the product cannot be disposed of as other household waste throughout the EU. To prevent potential damage to the environment or human health from uncontrolled waste disposal. Recycle responsibly to promote the sustainable use of material resources. If you want to return a used device, use the drop-off and collection system, or contact the retailer from whom you bought the product. The retailer can accept the product for environmentally safe recycling.

A declaration by the manufacturer that the product complies with the requirements of the applicable EU Directives.

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