
EN: SMILY – USER MANUAL

Product structure

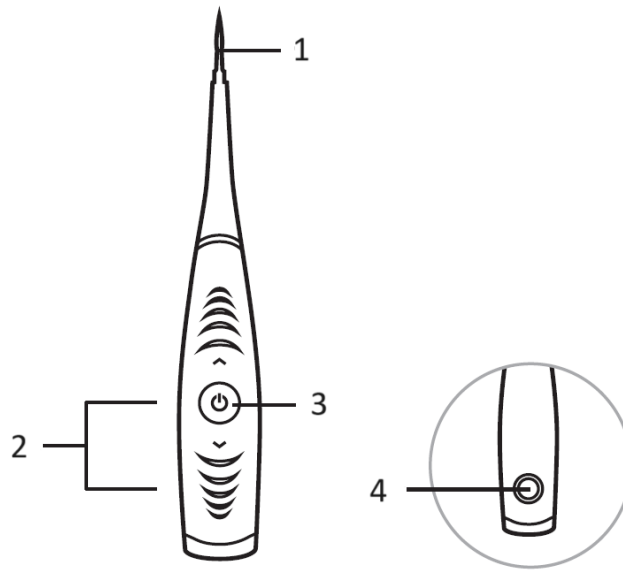1. Cleaner Head
2. +Key/-Key
3. ON/OFF button
4. DC Charger

This product is small in shape and convenient to carry. It allows you to clean your teeth at home. This product uses a clean head with the same structure as dentist. It can effectively remove obvious dental calculus, and also have certain cleaning effect on dental plaque and tooth stains.

Product features

- Professional grade alloy metal cleaning head, cleaning effect very nice;
- Five adjustable, suitable for different needs;
- Intelligent high frequency vibration, effective removal of dental calculus and so on;
- Whole body design, beautiful and generous;
- The battery is durable and practical.

Safety warnings

- Metal head of the tooth cleaner, not replaceable
- Do not place products in direct sunlight. Put at -10 C to 40 C indoor environment
- Do not wash or soak the product for a long time
- Products are restricted to adults only. Please keep them property after use, avoid children's contact, so that children will not be injured by mistake.
- The product contains batteries. Do not place products in a fire or high temperature container to avoid fire or explosion
- In case of product quality problem, please do not disassemble yourself. Please return to original plant for repair


Product Specifications

- Rated voltage : 3.7V
- Rated current : 130mA
- Battery capacity : 300mAh
- Standby current : 520uA
- Product size : 183.6x30x22.3mm
- Product weight : 53.6g
- Charging voltage : 5V
- Motor noise : 560dB

Usage

Before first use, the product must be charged fully (full charge: 2 hours). For the location of the charging port, please refer to number 4 under Product structure.

Press the boot button to turn on, adjust the key position through the upper and lower adjustment buttons, and start cleaning the parts that need cleaning (reference picture below). After cleaning the teeth, wash the metal head with flowing clean water. Clean the Department to wipe the body of the tooth cleaner.

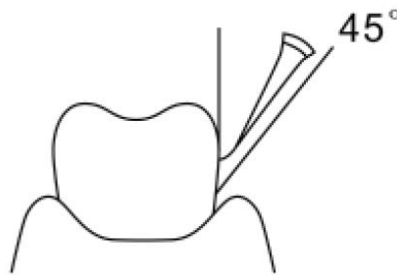Charger

1. This product uses the DC interface charging mode, the user takes out the distribution DC interface charging line, the DC interface is inserted into the machine back socket, and the USB interface is inserted at the end of the charger.
2. After entering the charging mode, the tooth cleaner enters the closed state. At this time, the machine will stop working
3. In charging, the light at the bottom is red and full of is green light (reference picture below).

INSTRUCTIONS FOR RECYCLING AND DISPOSAL:


This label means that the product cannot be disposed of as other household waste throughout the EU. To prevent potential damage to the environment or human health from uncontrolled waste disposal. Recycle responsibly to promote the sustainable use of material resources. If you want to return a used device, use the drop-off and collection system, or contact the retailer from whom you bought the product. The retailer can accept the product for environmentally safe recycling.

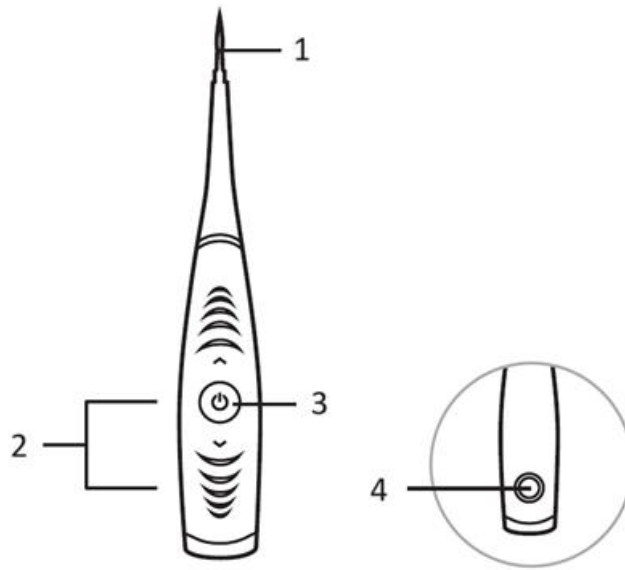A declaration by the manufacturer that the product complies with the requirements of the applicable EU Directives.

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